

## SecureSeal™ THERMAL ADHESIVE FILMS AND FOIL

### T329-1 & -2

#### SecureSeal™ THERMAL ADHESIVE SEALING FILM

for PCR application


This transparent sealing tape consists of a 2.0 mil polyolefin film coated on one side with a pressure sensitive acrylate adhesive which does not interfere with cycle reactions. It is ideal for reducing well-to-well contamination and/or spill over in sensitive PCR applications where the minimization of evaporation and vapor loss is critical.

SecureSeal™ Thermal Film was developed with the assistance of a major cycler manufacturer for PCR applications. Not only does it offer low-autofluorescence but it will prevent vapor loss and is thermostable and functional from -70 °C to 100 °C. Manufactured RNase, DNase and DNA-free. DMSO resistant.

Packed in a bag of hundred sheets.

**Note: Performance may depend upon the specific collection/sample vessel used as well as the specific conditions subjected to.**

Cat. #	Color	Sterile	Qty/Cs
T329-1	Transparent	No	1000
T329-2	Transparent	Yes	1000

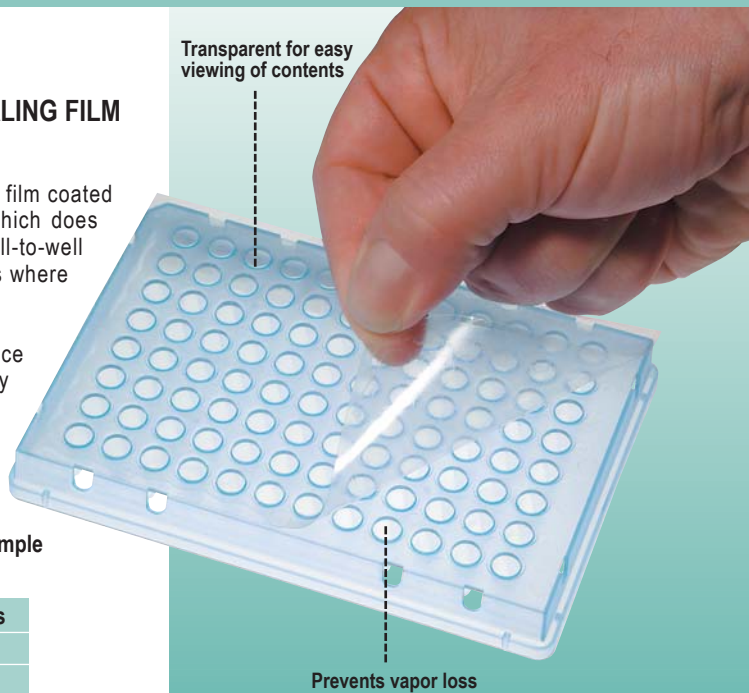

### T329-5

#### SecureSeal™ ALUMINUM SEALING FOIL

This type of material is ideal for manual sealing during PCR work and also for high throughput applications. Adhesive backing makes it easy to apply. Will resist temperatures from -86 °C to 110 °C. It is recommended to use the Amplate™ Roller (Cat. no.T329-9) for ensuring a perfect bond, eliminating the danger of evaporation. Pierceable with a pipet tip for easy access to sample. DMSO resistant.

Cat. #	Description	Qty/Pk
T329-5	Peeling foil	100 sheets



### T329-9

#### AMPLATE™ ROLLER

For ensuring a perfect seal when using either SecureSeal™ Thermal sealing film or foil on PCR plates. Roller made of medium hard rubber. Heavy-duty handle with comfort grip reducing fatigue. Will last a long time.

Cat. #	Size	Qty/Pk
T329-9	10.16 cm (4 in.)	1



### T329-10

#### AMPLATE™ MAT

Made of thermoplastic rubber

This flexible sealing cover is used on 96-well plates along with clip down and screw top thermal cyclers and has been proven to be a secure and effective way of sealing. Since it is reusable, it is a nice way to make this step of the procedure cost effective. Dimples on one side of the mat ensure it is well placed over the tubes. Autoclavable.

Cat. #	Color	Qty/Pk
T329-10	Blue	5